


FACIAL FAT TRANSFER SET

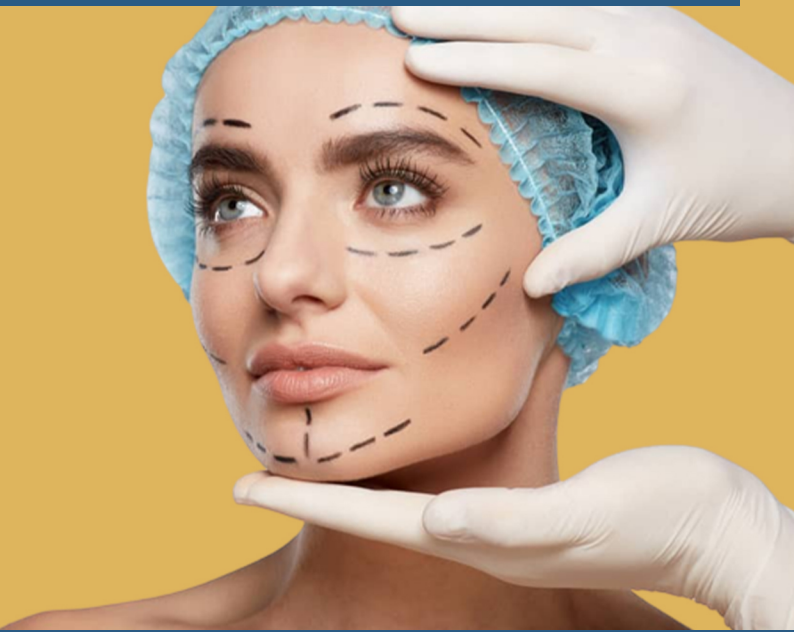

Redefining the Art of Fat Grafting

Facial Fat Transfer Set

Alira's Facial Fat Transfer Set is an extensive collection of reusable and cell-friendly cannulas for redefined facial sculpting. Designed by recommendation of experts, the set is equipped with all instruments to achieve a natural and harmonious result when sculpting and contouring the face, rejuvenating the skin and revising scars.

1 Tumescant Infiltrator Straight

The safe and effective administration of tumescent solutions using tumescent infiltration cannulas ensures a painless liposuction procedure.



Length: 20cm | **Diameter:** 2.5mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub | Power Cannula | Fixed Handle | PAL Handle

2 Aggressive Multiport Harvester Cannula - One End Sharp

Alira's Aggressive Harvester Cannula – One End Sharp is essential for liposuction procedures or fat harvesting. A harvester designed for speed. Fully equipped with beveled and elevated sharp ports positioned asymmetrically to create a helical flux.



Length: 15cm | **Diameter:** 2.5mm | **Qty:** 01

Length: 20cm | **Diameter:** 2.5mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub | Power Cannula | Fixed Handle | PAL Handle



3 Carraway Style Cannula - Three Spiral Port

Alira's Carraway Harvester is a spiral three-port design that is popular for small volume treatments carried out under local anesthesia and is safe and efficient for multi-plane adipose tissue harvesting.



Length: 15cm | **Diameter:** 2mm | **Qty:** 01

Length: 20cm | **Diameter:** 2.5mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub
| Power Cannula | Fixed Handle | PAL Handle

4 Blunt Tip Injector Cannula - Straight

Alira's featured Blunt-Tip Injector Cannula passes through soft tissue effortlessly. In addition, this cannula is available in different sizes and working lengths depending on surgical preference. You can also choose between Luer-Lock Type connector, Super Luer Lock Type Connector or Threaded Type connector.



Length: 5cm | **Diameter:** 0.7mm | **Qty:** 01

Length: 7cm | **Diameter:** 1.2mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub
| Power Cannula | Fixed Handle | PAL Handle

5 V Shaped Dissector Cannula-Single Port

This cannula is specifically designed for use in the superficial compartment. The cutting notch is very useful for releasing dense fibrous tissue. Such as acne scars and nasolabial folds. In addition, the cannula is available with different sizes and working lengths depending on surgical preference. You can also choose between Luer-Lock Type connector, Super Luer Lock Type Connector, Threaded Type connector and with Fixed Power Handle as per your surgical needs.



Length: 10cm | **Diameter:** 1.5mm | **Qty:** 01

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub
| Power Cannula | Fixed Handle | PAL Handle

6 Little Injector Cannula

In dense and fibrous tissue, a blunt tip of Alira's Little Injector cannula with an open, pointed tip works very well.



Length: 4cm | **Diameter:** 0.7mm | **Qty:** 02

Length: 5cm | **Diameter:** 0.7mm | **Qty:** 02

Length: 5cm | **Diameter:** 0.9mm | **Qty:** 02

Connector Type Options: Luer Lock | Ultra Luer Lock | Nouvag Threaded Hub | 60cc Toomey Hub
| Power Cannula | Fixed Handle | PAL Handle

7 Snap Locks

In order to provide negative syringe pressure during liposuction procedures, Alira's Snap Lock is a special tool. The liposuction syringe's plunger flanges may accommodate these syringe Snap Locks.



Type: 20cc | **Qty:** 01

Type: 60cc | **Qty:** 01

8 Johnnie Lock

Johnnie Lock holds the plunger under vacuum pressure for easy syringe harvesting. The plunger is drawn back to the desired pressure level then twisted and locked securely.



Type: 10cc | **Qty:** 01

9 Anaerobic Transfers

Alira's Anaerobic Transfers utilizing Micro Injector Cannulas are used to prepare and transport harvested tissue for fat grafting.




Type: Luer to Luer 2.4mm | Qty: 02

Type: Luer to Luer 1.4mm | Qty: 01

Type: Luer to Luer 1.2mm | Qty: 01

10 Decant Stands

To prepare and transfer harvested tissue for reinjection, anaerobic transfers are performed. Alira's Decant Stand's design also enables the decant of extracted tissue using natural gravity.

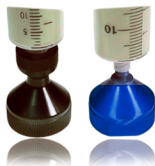

Type: 60cc to Luer | Qty: 01

Type: Luer Lock to Luer Lock | Qty: 01

11 Mesh For Nano Transfer

For separating big adipose tissue from small re-graftable fat cells through a Nano filter, Nano fat transfer mesh is excellent.



Size: 300 Micron | Qty: 01

Size: 500 Micron | Qty: 01

Size: 800 Micron | Qty: 01

12 Nano Transfer

Alira's Nano Transfer is a tool made to uniformly enlarge collected adipose tissue so that microneedles can readily inject it. The Nano Transfer has a cartridge. When doing liposuction and fat transfer operations, it is utilized to filter fat.



Qty: 01



WHY CHOOSE ALIRA?

Enhanced Precision: You can achieve precise fat graft placement, resulting in consistent and desirable cosmetic outcomes.

Minimized Patient Trauma: Alira's instruments minimize tissue trauma, leading to faster healing and reduced post-operative discomfort for patients.

Reduced Surgeon Fatigue: Ergonomic handles improve your comfort during procedures, reducing the risk of fatigue and enhancing surgical control

Biocompatible Materials: All materials used are biocompatible, ensuring compatibility with the human body and reducing the risk of adverse reactions.



CORE AREAS OF EXPERTISE

Manufacturing
Extensive Range of
High Quality
Surgical Instruments

Instruments
Customization

Instruments
Repairing Services

Contract
Manufacturing

CONTACT US

ALIRA MEDICAL DEVICES LLC

📞 +1 (818) 614 4592 ✉ info@alirameditdevices.com 🌐 www.alirameditdevices.com

📍 14622 Ventura Blvd Ste 102, #712 Sherman Oaks, CA 91403

